

Resonator Design and Capacity

Zach Martineau, PE
Senior Mechanical Engineer




Technical
InterChange 2026

1

Resonator Design and Capacity



- Resonator Design

- Resonance
- Three Masses

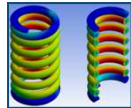

- Resonator Capacity

- Classifying Payload Mass
- LabRAM II Mix Capacity

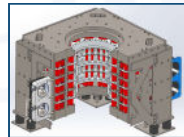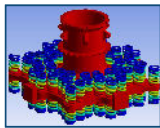


2

Technical
InterChange 2026

2

Resonator Design- Power Obstacle

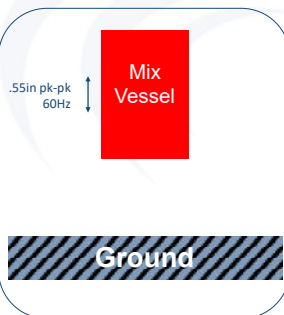

- Oscillating Motion Causes Mixing

- ≈60 Hz
- 100g of acceleration

- How to Power It?

- 5 lbm – 2,600 W (3.5 HP)
- 375 lbm – 195 KW (260 HP)
- 3000 lbm – 1.56 MW (2150 HP)




3

Technical
InterChange 2026

3

Resonator Design- Payload in Resonance

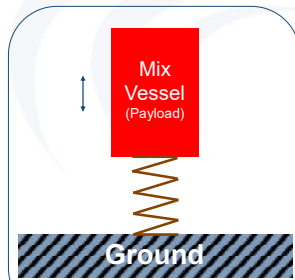

- **Solution**
 - Add a spring
 - Resonant Condition:



$$f \propto \sqrt{\frac{k}{m}}$$

- **Problem**
 - Astronomical force to ground
 - 375 lbm payload results in $\pm 37,500$ lbf (167kN)


Resodyn
ACOUSTIC MIXERS

4

Technical
InterChange 2026

4

Resonator Design- Cancellation Mass

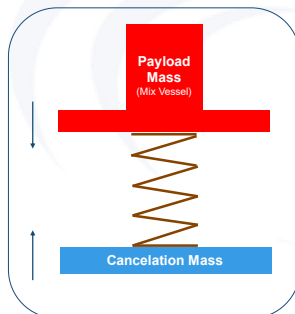

- **Solution**
 - Add a cancellation mass

- **Problem**
 - How do you hold on to it?


Resodyn
ACOUSTIC MIXERS

5

Technical
InterChange 2026

5

Resonator Design- Cancellation Mass

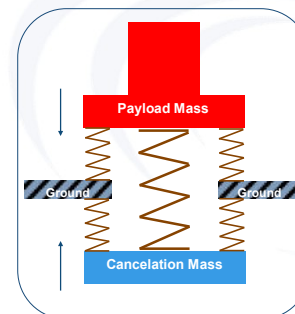

- **Solution**
 - Add springs to ground between the masses

- **Problem**
 - Payload variation means the system isn't always perfectly tuned resulting in force to ground


Resodyn
ACOUSTIC MIXERS

6

Technical
InterChange 2026

6

Resonator Design- Frame Mass

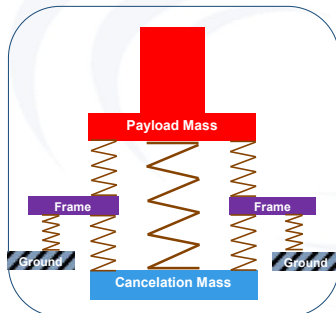

- **Solution**
 - Add a third mass (frame) to cancel the remaining force to ground
 - Maximum alternating force to ground from a RAM55 is ± 80 lbf (356N)


Resodyn
ACOUSTIC MIXERS

7

Technical
InterChange 2026

7

Resonator Design

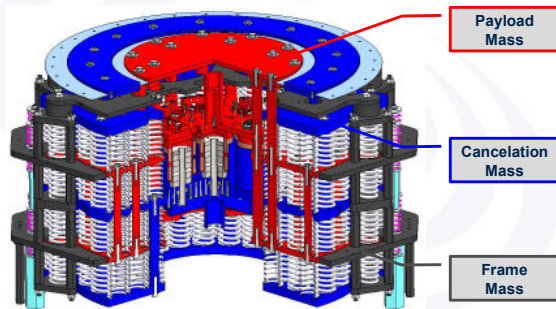

Resodyn
ACOUSTIC MIXERS

8

Technical
InterChange 2026

8

Mixer Capacity



Model	Capacity
LabRAM I	1.1 lbm / 0.50 kg
LabRAM II	2.2 lbm / 1.0 kg
OmniRAM	11 lbm / 5.0 kg
RAM 5	80 lbm / 36 kg
RAM 55	924 lbm / 420 kg

Resodyn
ACOUSTIC MIXERS

9

Technical
InterChange 2026

9

Classifying Payload Mass: Dead Mass

- **Dead Mass: Mass That Does Not Absorb Energy**

- Vessel Holder
- Vessel
- RTD Probes
- Fasteners
- Components that do not move relative to the payload plate



Resodyn
ACOUSTIC MIXERS

10


Technical
InterChange 2026

10

Classifying Payload Mass: Live Mass

- **Live Mass- Mass That Absorbs Energy**

- Powder
- Liquids
- Paste
- Coolant
- Mass that moves relative to the payload plate

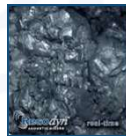

Resodyn
ACOUSTIC MIXERS

11

Technical
InterChange 2026

11

Classifying Payload Mass: Total Equivalent Mass

$$M_{Eqv} = M_{Dead} + M_{Live} \times C$$

- C- Coupling Coefficient (Determined Experimentally)
- Coupling Factors:
 - Mix material properties
 - Vacuum
 - Headspace
 - Vessel geometry
 - And more...

Resodyn
ACOUSTIC MIXERS

12

Technical
InterChange 2026

12

Classifying Payload Mass: 100g Equivalent Mass Limits

Model	Capacity	Maximum 100g Equivalent Mass
LabRAM I	1.1 lbm / 0.50 kg	4.2 lbm / 1.9 kg
LabRAM II	2.2 lbm / 1.0 kg	5.5 lbm / 2.5 kg
OmniRAM	11 lbm / 5.0 kg	80 lbm / 36 kg
RAM 5	80 lbm / 36 kg	375 lbm / 170 kg
RAM 55	924 lbm / 420 kg	3100 lbm / 1400 kg

$M_{Eqv} = M_{Dead} + M_{Live} \times C$



13



13

Classifying Payload Mass: Examples



Model	Payload Classification	Maximum 100g Equivalent Mass
LabRAM II	2.2 lbm	5.5 lbm

	Standard Vessel Holder 8oz PP Jar	2.4 lbm
	Oversized Vessel Holder 32oz Glass Jar	6.8 lbm
	Standard Vessel Holder 8oz PP Jar 1.0 lbm Sand	3.4 lbm




14



14

Classifying Payload Mass: Examples

- Increasing the Dead Mass Lowers the Operating Frequency but does not Change the Mix Power

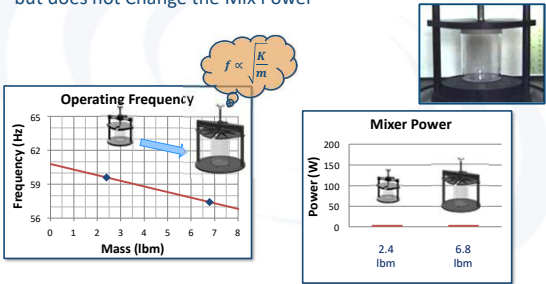


15



15

Classifying Payload Mass: Examples

- Adding Live Mass Changes both the Operating Frequency and Mix Power

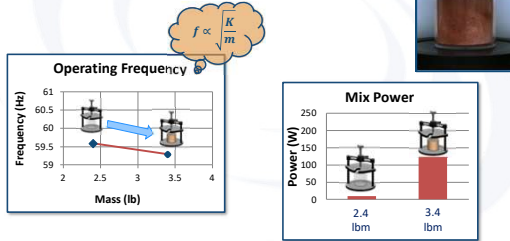


Resodyn
ACOUSTIC MIXERS

16

Technical
InterChange 2026

16

LabRAM II Capacity: Operating Conditions

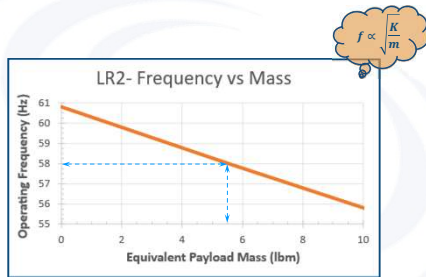

$$M_{Eqv} = M_{Dead} + M_{Live} \times C$$


Resodyn
ACOUSTIC MIXERS

17

Technical
InterChange 2026

17

LabRAM II Capacity: Acceleration De-Rate Conditions

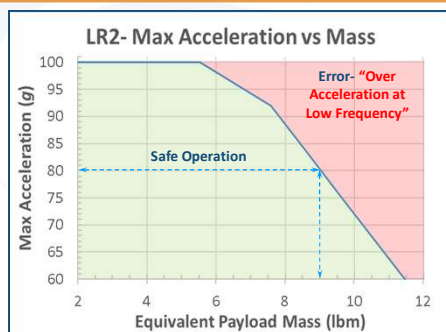

Resodyn
ACOUSTIC MIXERS


18

Technical
InterChange 2026

18

LabRAM II Capacity: Available Power

- Work = Force x Displacement
- Usable Power at 100g: 600 W
- For Example:
 - If a live mass absorbs $250 \frac{W}{lbm}$
 - $\frac{600 W}{250 \frac{W}{lbm}} = 2.4 lbm \text{ Max}$
 - Mixing more will result in error-
"Amplifier Over Current"

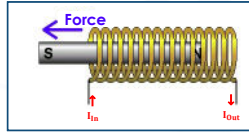

19

Mixer Capacity: Available Power



Model	Maximum Available Power
LabRAM I	230 W
LabRAM II	600 W
OmniRAM	4.2 kW
RAM 5	24.7 kW
RAM 55	235 kW

20

Mixer Design and Capacity- Summary

- Resonance is Paramount to Efficient RAM Operation
- Three Masses are Required to Manage Floor Vibrations
- Payload Mass (Equivalent Mass) is:

$$M_{Eqv} = M_{Dead} + M_{Live} \times C$$

- Resonators Have Both Mass and Power Limits

21

Thank you for your time
and attention.





22

Technical

InterChange 2026

22